HANDPIECES

July - August 2025

















🙆 in

f

(Dublin 01 427 3700/Cork 021 602 0544

www.dmi.ie

info@dmi.ie



Dublin: 01 427 3700 /Cork: 021 602 0544

0544 Web: www.dmi.ie

Email: info@dmi.ie

All prices exclude VAT. E&OE.

f) 🞯 in

Dentsply Sirona

5

47% OFF





Synea Vision WG-56LT

Optic

- 1:1
- 1 Port Spray
- Stainless Steel 2 Year Warranty

RRP €1,377 €708

WH-30009000



 Stainless Steel 1 year warranty RRP €424



€165

CA 1:1L Bien Air Dental







Dublin: 01 427 3700 /Cork: 021 602 0544

Email: info@dmi.ie

All prices exclude VAT. E&OE.



Sirona Tl Energo 1:5 L

- Dentsply Sirona Longer service life thanks to Diamond like Carbon coated ball bearings
- Optimal illumination for an excellent view
- 1:5 ratio has four hole spray
- 3 year warranty

SIR-6812056

SIR-6812148 Sirona T1 Energo S 1:5 L - SHORT

RRP €1.381 Now €737



Sirona T2 Energo 1:5 L

- Robust ceramic ball bearings
- Pleasant and easy to grip thanks to its titanium sleeves
- Optimal illumination for great
- visibility during treatments
- 1:5 ratio has four hole spray
- 2 year warranty

SIR-6812080

RRP €922 Now €492

Sirona T3 Energo Dentsply Sirona 1:5 L - NON OPTIC Contra-Angle Red • Ratio 1:5 • Non optic • 1 year warranty SIR-6812114 RRP €697 47% OFF Now €372

EXPERTmatic E25L



KAVO

MASTERmatic M25L

- Optic
- Ratio: 1:5
- 3-Port Spray
- Plasmatec Anti-Slip Finish
- Perfect grip and excellent hygiene
- 2 Year Warranty

KAVO-M25I RRP €2,198

Only €1,260 each

MASTER Handpiece Range Mix & match - Triple Pack only €833 each

Ti-Max Z95L Switch

- Optic contra-angle
- 1:5 Increasing
- Quattro spray
- Ceramic ball-bearings World first switch between
- water and spray on handpiece 3 year Warranty

NSK-C1034001SW RRP €1.866

€1,493 eachor Buy 3 @ only €995 each **Plus FREE**

NSK Z Micro (Ultra Mini Head) Handpiece FREE worth €1,346 (4 handpieces in total)

NEW Synea Fusion

Power

Designed specifically for removal of zirconia and other hard materials

- Optic
- 1:4 Speed Increasing
- 3-Port Spray
- Scratch-Blocker coating
- 3 year warranty

WH-30459000



Optic

BIEN-1600690-001 Micro Series BIEN-1600386-001 Standard

Web: www.dmi.ie

Email: info@dmi.ie




Dental

BIEN-1600940-001 Micro Series BIEN-1600941-001 Standard

Buy 4 @ only €600 each




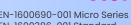







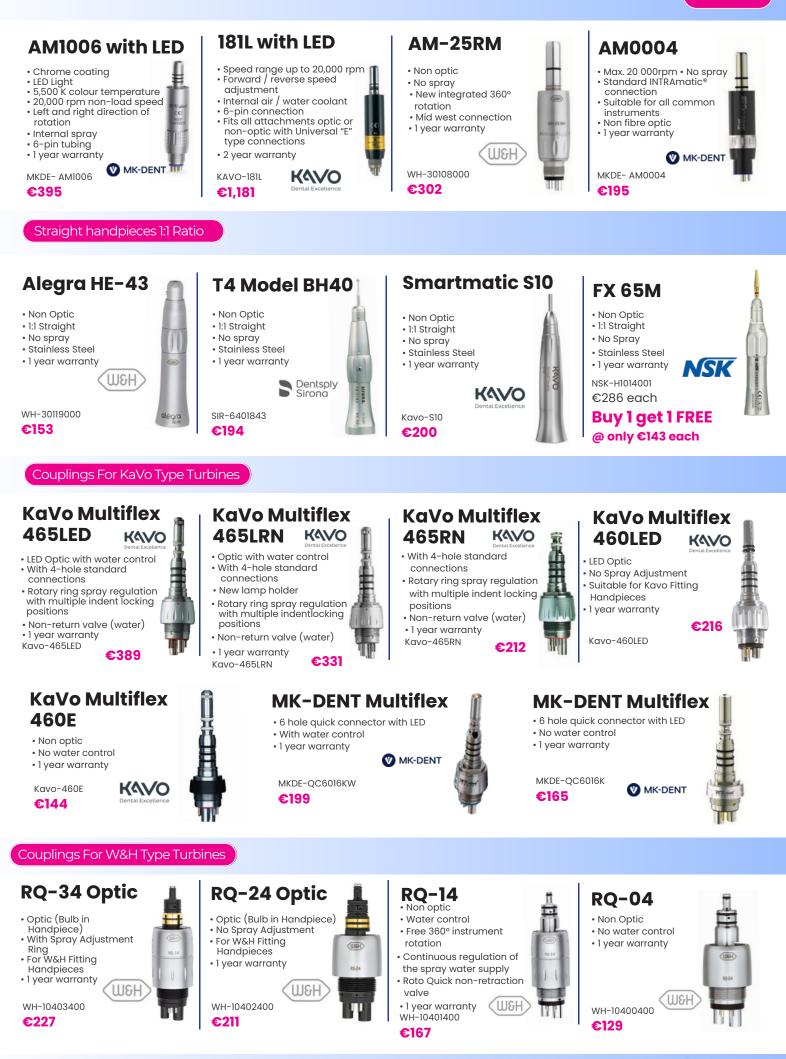

Dublin: 01 427 3700 /Cork: 021 602 0544

Web: www.dmi.ie

Email: info@dmi.ie

All prices exclude VAT. E&OE.

f) 🞯 in



Email: info@dmi.ie

All prices exclude VAT. E&OE. 👔 👩 ท

Ti-Max X-SG65



- Non Optic
- Gear Ratio : 1:1 Direct Drive
- For HP Burs (Ø2.35)
- Water Spray : External Cooling
- Body Material : Titanium
- Body Coating : DURACOAT. Max Speed : 40,000 min-1 NSK-H1038001

Ti-Max X-SG65L



- Gear Ratio : 1:1 Direct Drive For HP Burs (Ø2.35)
- Water Spray : External Cooling.
- Body Material : Titanium
- Body Coating : DURACOAT
- Max Speed : 40,000 min-1

NSK-H1009001





- Non Optic
- Gear Ratio : 20:1 Reduction. Water Spray : External and Internal Cooling (Kirschner and Mever)
- Body Material : Stainless Max Speed : 2,000 min-1
- Max Torque : 80 Ncm


NSK-C1010001

Ti-Max X-SG20L



- Optic
- Gear Ratio : 20:1 Reduction Water Spray : External and Internal
- Cooling (Kirschner and Meyer) Body Material : Titanium Body Coating : DURACOAT. Max
- Speed : 2,000 min-1 Max Torque : 80 Ncm
 - NSK-C1003001

Implantology & Oral Surgery

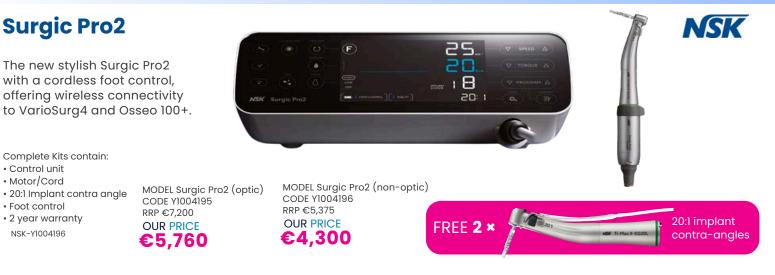

MYLUNOS DUO

MyLunos Duo combines two of the most important treatment methods for a prophylaxis session: powder jet treatment and ultrasound treatment. You can use the functions Supra and Perio to switch as required between supragingival and subgingival powder jet application. Here, two different Lunos® Prophy Powders can be easily interchanged.

Exclusive to DMI -Plus Free User Training

Clearly laid out and

hygienic display

٠



MYLUNOS Pro Powder Jet Handpiece

For supragingival and subgingival applications


Contents:

- MyLunos® powder jet unit ٠ including Supra nozzle
- Perio nozzle
- 8 Perio tips
- 2 powder containers with tray
- Lunos® Prophy Powder Gentle Clean and Perio Combi
- O-ring set

RRP. €2,395 OUR PRICE £1,

DURR-2034910053 Bien Air DURR-2034910050 KaVo DURR-2034910054 NSK



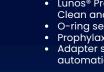



Woodpecker O - Star The Ultimate Curing Expert

- Performs light curing in 1 second
- Longer wavelength: 385 515 nm
- It includes a light intensity meter which is located in the charging base
- It has 7 different work modes: normal, high intensity, turbo, ortho, soft, pulse and verification
- Light intensity: up to 3000mw/cm²
- Small, rotating head allows you to reach all angles of the mouth and reach the back teeth without difficulty

WOOD-OSTAR

€495

Woodpecker LED.H ORTHO

Specially designed for Orthodontics

- Specially designed for orthodontics. In just 3 seconds the brace solidifies
- Short, powerful cycles without overheating
 Constant work without losing intensity, for a
- complete curing of the resinsLight and easy to use, with a rechargeable
- battery
 The new bevelled light guide design can easily reach any tooth and radiate vertically at any angle
- Includes a light intensity meter

WOOD-LED.H.ORTHO



WOODPECKER

Light

WOODPECKER

Woodpecker iLED

Curing Light

- Ergonomic curing light with 360° rotating head
- Fast Cures up to 2mm in just one second
- Light intensity up to 2300mw/cm² with constant power output
- Power input: AC100V-240VLight output intensity: 1000
- mW/cm2~2500mW/cm2 • Time Sets: 5-20 Seconds
- Time sets: 5-20 seconds
 Operation Mode: turbo mode, normal mode





WOODPECKER

WOODPECKER

Woodpecker LED.H

• Woodpecker LED-H Light Curing Lamp provides reliable light curing of any type of material

€295

- It is a light weight lamp, which will allow you to use it during a large part of your professional daily activity, without causing fatigue
- Three working modes: Normal, ramp or intermittent
- High-capacity battery with charge duration for 200 continuous lightcuring cycles

WOOD-LED.H



Shake the can 3 or 4 times, hold the

can in upright position, then gently

press in the push button to directly

spray into the bur insertion hole.

5

Handpiece Maintenance Tips

How to Properly Oil the Chuck of Dental Handpieces

The chuck is a critical part of your dental handpiece—it holds the bur securely in place and endures high-speed rotation and heat. Proper lubrication ensures optimal performance and longevity.

Why Chuck Lubrication Matters:

- Prevents corrosion and sticking
- Maintains grip strength on burs
- Reduces wear and failure of internal components

Remember to oil the chuck once a week for turbines and blue bands and before every autoclave cycle for red bands as these are really fast and need extra oil.

Why You Should Not Use Wipes—Especially Alcohol Based—on Dental Handpieces

While it might seem convenient, wiping your handpieces with disinfectant wipes-particularly those that are alcohol-based-can cause serious long-term damage.

What's the Risk?

Alcohol Breaks Down Seals and Gaskets:

- Alcohol and other harsh disinfectants can dry out and crack internal O-rings and seals, especially in air and water line connections
- This leads to air leaks, reduced performance, and eventual failure
- Degrades handpiece coatings

Invalidates Warranties:

· Many handpiece manufacturers specifically state that wiping with alcohol-based products will void the warranty

What to Do Instead:

- Use Manufacturer-Approved cleaning solutions
 - Follow the correct cleaning protocol:
 - Pre-clean with a damp gauze or cloth (water only or approved detergent).
 - Lubricate per the manufacturer's instructions
 - Sterilize in an autoclave







Scan the QR code to watch a video on the correct way to oil the chuck.